Conditioning

artners

CONDITIONING APK 1200C CANISTER CLEANER


CANISTER CLEANER

The Canister must be cleaned for re-use. If it is used again without cleaning, the remained substances will affect the result of analysis. So that, it must be cleaned by high purity gas on high temperature. This instrument can simultas neously condition 4 channels of canisters.



FEATURES

- Save conditioning time with 4 canister simultaneous cleaning
- Temperature setting for each channel
- Repeatedly pressurized and depressurized using high purity purge gas to condition
- Reduce internal pressure and enable immediately when completed conditioning.
- Cleaning methods such as vacuum, pressurization, temperature are easy to use with dedicated programs

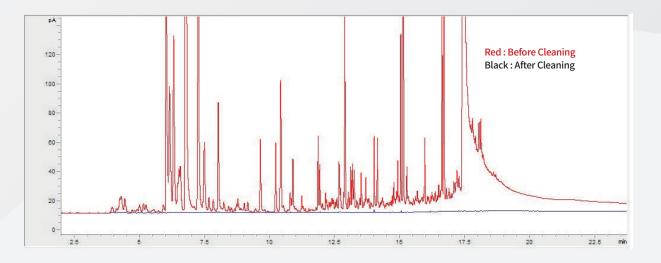

SPECIFICATION

SPECIFICATIONS	
Canister Heating	Room temp + 20°C~150°C(Max)
Cleaning Canisters	4 Canisters, 6L
Vacuum Pump	Rotary Type (5x10 ⁻³ Torr)
Tube Material	Flexible Stainless Steel Tube
Control	Touch screen control
Dimension	880 x 720 x 770 (W x D x H mm)
Power	220VAC, 50/60Hz, Max. 2.2KW, 10A